

Regd. Office :


A- 6, NAND Complex,
Near AIMS Oxygen, Juna Padra Road,
Baroda 390 020 - INDIA

Marketing / Customer Support :

5-6, 3rd Floor, Mrudula Sadan,
Pratap Road, Raopura,
Baroda 390 001 - INDIA

JUTRON Vision

enhancing vision ~ enriching life

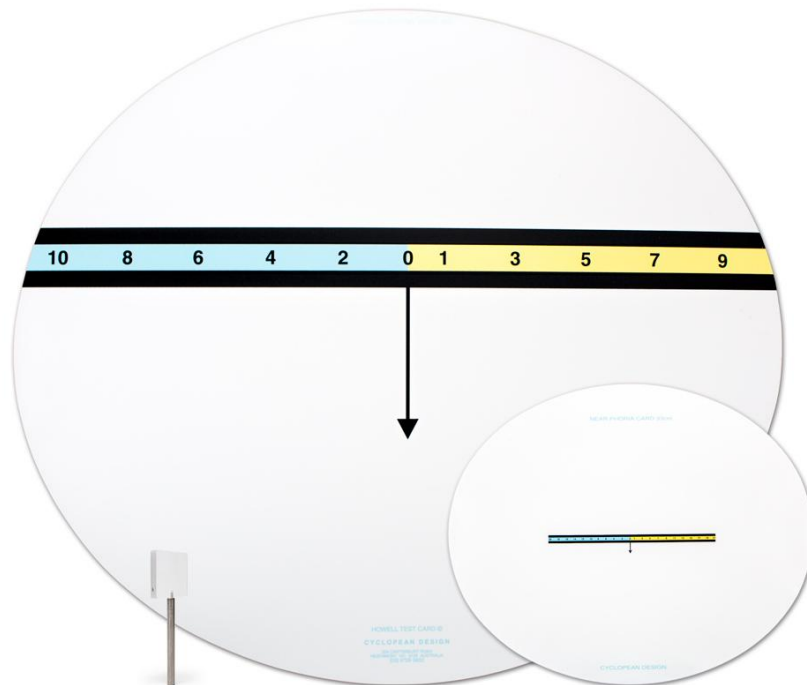

CDHP+ Howell Phoria Card™

Howell Phoria Card™

Both Near and Far come with a single 6 PD stick prism.

These cards allow direct measurement of the distance and near phorias in real space or with a phoropter. The patient simply has to tell which number the arrow points to. Any odd number (in yellow) shows esophoria while the even numbers (in blue) shows exophoria. An AC/A can quickly be determined with trial lenses.

The reverse side of the near card has a X-Cyl target, and a variety of letter and word targets.

Near card is about 12" wide.

Distance chart is about 24" wide.