

Regd. Office :


A- 6, NAND Complex,
Near AIMS Oxygen, Juna Padra Road,
Baroda 390 020 - INDIA

Marketing / Customer Support :

5-6, 3rd Floor, Mrudula Sadan,
Pratap Road, Raopura,
Baroda 390 001 - INDIA


JUTRON Vision

enhancing vision ~ enriching life

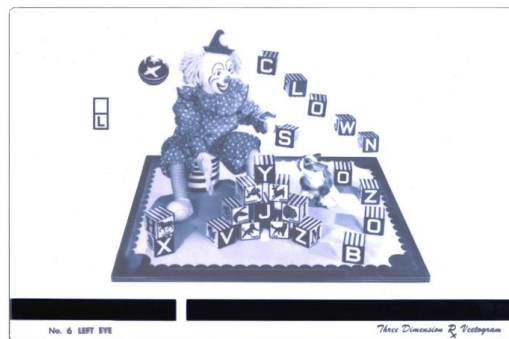

SOV6 Variable Vectogram (Clown Vectogram)

Variable Vectogram Clown Vectogram

- 8-1/4" x 5-1/2"
- Polarized vectograms ensure testing and training at actual distances. Vectograms stabilize fusion and stereopsis, eliminate suppression and develop simultaneous vision.
- Used in monocular, binocular pursuit and saccadic training, vectograms help correct and prevent anomalies involving projection and hand-eye coordination.
- Ideal for use with Bernell's Dual Polachrome™ Trainer
- Polarized goggles included

Cleaning Instructions:

DO NOT IMMERSE THE VECTOGRAM PANELS IN WATER

To clean the vectogram panels and guide, use a soft, damp, lint-free cloth. Dampen cloth using glass cleaner or mild detergent/water.

If panels are removed from guide during cleaning, replace the panels in the guide, placing the panel with the solid bar on top of the panel with the numbers/letters and ensure that the plastic portion of the guide is behind the panels

Storage Instructions:

Store Vectogram in a dry place, away from direct sunlight.

Warranty: 1 year limited manufacturer warranty from the date of purchase

Delhi : Chandigarh : Lucknow : Jaipur : Indore : Ahmedabad : Mumbai : Pune : Nagpur : Hyderabad : Bangalore : Chennai

Tel. : +91-99987-02020 • mail@jutronvision.com • www.JutronLowVision.com • www.JutronVision.com